
SINCE 1966
MADE IN ITALY

Air Pro

by **FLAEM**® **3000**
Plus Version

- I** ISTRUZIONI PER L'USO - Aerosol professionale per uso ospedaliero, home care e noleggio
- GB** INSTRUCTIONS FOR USE MANUAL - Professional aerosol therapy unit for hospital use, home care and renting
- F** MODE D'EMPLOI - Générateur d'aérosol professionnel à usage hospitalier, de soin à domicile et de location
- NL** GEBRUIKSAANWIJZING - Professionele aerosol voor gebruik in ziekenhuizen, home care en verhuur
- D** BEDIENUNGSANLEITUNG - Professioneller Inhalator für den Gebrauch im Krankenhaus, Homecare und zum Mieten
- E** MANUAL DE INSTRUCCIONES - Nebulizador profesional para uso hospitalario, asistencia domiciliaria y alquiler
- RUS** ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ - Ингалятор для профессионального использования в больницах, для использования в домашних условиях, а также для предоставления напрокат
- GR** ΕΓΧΕΙΡΙΔΙΟ ΟΔΗΓΙΩΝ ΧΡΗΣΗΣ - Επαγγελματικό αερόλυμα για νοσοκομειακή χρήση, φροντίδα στο σπίτι και ενοίκιαση

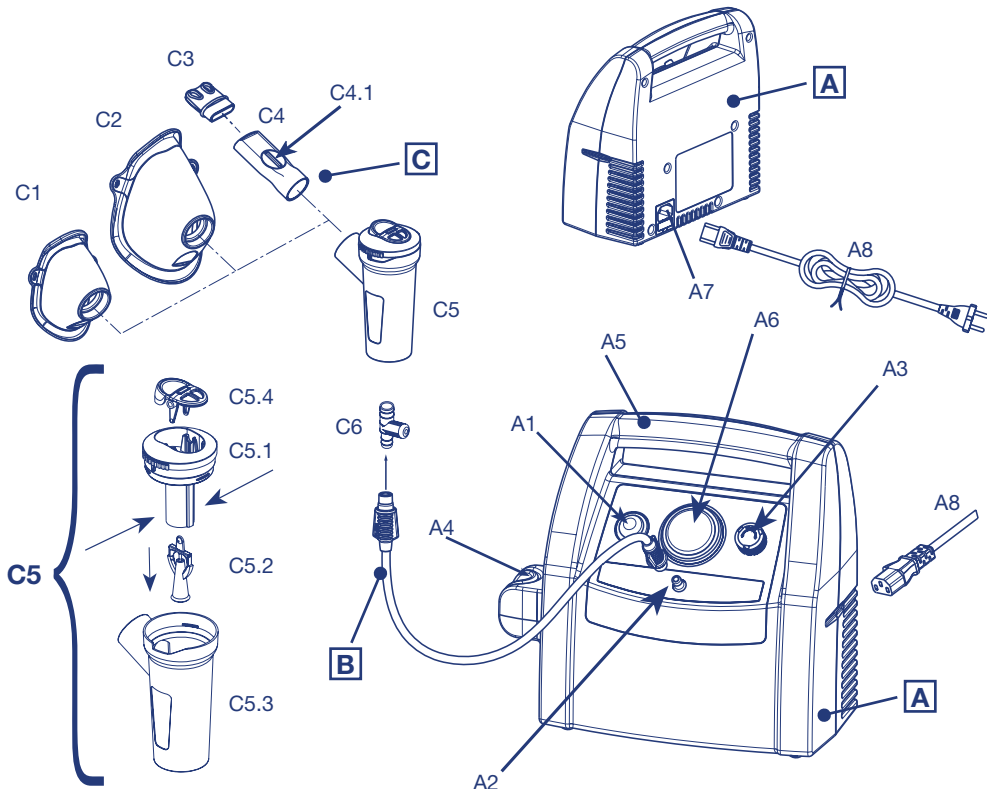SCHEMA DI COLLEGAMENTO PER NOMENCLATURA - CONNECTION DIAGRAM FOR NOMENCLATURE
 SCHEMA DE CONNEXION POUR NOMENCLATURE - NOMENCLATUUR VERBINDINGSSCHEMA
 SCHALTPLAN FÜR DIE BEZEICHNUNGEN - ESQUEMA DE CONEXIÓN PARA NOMENCLATURA
 СХЕМА СОЕДИНЕНИЯ СОСТАВНЫХ ЧАСТЕЙ - SCHEMAT POŁĄCZEŃ DLA NOMENKLATURY
 ΣΧΕΔΙΑ ΣΥΝΔΕΣΗΣ ΓΙΑ ΟΝΟΜΑΤΟΛΟΓΙΑ

Air Pro 3000

Mod. P0712EM F3000 Plus Version

Apparecchio per aerosolterapia

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Mod. P0712EM F3000 Plus Version

Aerosol therapy apparatus

We are pleased you have purchased our product and we thank you for your trust in us.

Read these instructions carefully and retain them for future reference.

Only use the apparatus as described in this manual.

ENGLISH

INTENDED USE

Air Pro 3000 Plus Version is a professional medical device whose intended use is the administration of medication prescribed or recommended by your doctor for the treatment of breathing ailments. It can be used in health care structures like hospitals, in home care and for hire. The product should be used by legally authorized medical personnel (doctors, nurses and therapists), by the home therapist, or by the patients themselves under the guidance of qualified medical personnel. It is highly recommended that the health professional and/or the patient read and understand the use and maintenance information of this manual.

Carefully follow the user and maintenance instructions to ensure proper functioning and to extend the life of the device.

Air Pro 3000 Plus Version is fitted with a knob to adjust the pressure (A3), and with pressure gauge (A6) to display the set pressure value.

The device does not require lubrication, is manageable, simple to use, reliable, resistant and quiet.

Note: Use only original FLAEM accessories

DEVICE NOMENCLATURE

A Aerosol Device (main unit)

- A1- Switch
- A2 - Air outlet port
- A3- Pressure adjustment knob
- A4- Nebulizer port
- A5- Carrying handle
- A6- Pressure gauge
- A7- Power supply cord socket
- A8- Power cord

C Accessories

- C1- Child mask
- C2- Adult mask
- C3- Non-invasive adult nasal prong
- C4- Mouthpiece with valve
 - C4.1- Expiratory valve
- C5- RF7 Dual Speed Plus Nebulizer
 - C5.1- Upper part
 - C5.2- Nozzle
 - C5.3- Lower part
 - C5.4- Speed selector with valve system
- C6- Nebulization manual control
- C7- Replacement valves for nebuliser
- C8- Air filter replacement kit (2pcs)

B Connection tube (main unit / nebuliser)

⚠ IMPORTANT NOTICES

The manufacturer makes every effort to ensure that every product is of the highest quality and safety; however, as for any electrical device, basic safety regulations must always be followed in order to avoid harming persons and things.

- Before using the product for the first time, and periodically during its lifetime, check the power cable to make sure it is not damaged. In the event it should be damaged, do not plug in the cable and immediately bring the product to an authorised FLAEM service centre or to your trusted dealer.
- Average estimated life of accessories is one year. However, in order to assure continued maximum therapeutic effectiveness, single-patient use is advisable and the nebuliser should be replaced every six months during intensive utilisation periods (or earlier if the nebuliser is clogged).
- In the presence of children and non- self-sufficient individuals, the device must be used under the close supervision of an adult who has read this manual.
- Some parts of the device are small enough to be swallowed by children; keep the device out of the reach of children.
- Do not use the provided tubing and cables for any other purpose than those specified, as they can cause risk of strangulation. Be particularly attentive with children and individuals with disabilities because they are often unable to correctly assess risk.
- The apparatus is unsuitable for use in presence of flammable anaesthetic mixture with air, oxygen or nitrous oxide.
- Always keep the power cable away from hot surfaces.
- Keep the power cable away from animals (for example, rodents) which could damage the insulation.
- Do not handle the plug with wet hands. Do not use the device in damp conditions (for example when taking a bath or shower). Do not immerse the device in water; in the event of immersion immediately disconnect the plug. Do not remove or touch the immersed device; unplug the power cable first. Immediately bring the device to an authorised FLAEM service centre or to your trusted dealer.
- Use the device only in dust-free conditions, otherwise treatment could be compromised.


- The device casing is not waterproof. Do not wash the device under running water or by immersion and keep it safe from being sprayed by water or other liquids.
 - Do not expose the device to particularly extreme temperatures.
 - Keep the device away from sources of heat, direct sunlight or excessively warm locations.
 - Do not obstruct or introduce objects in the filter and its relevant housing in the device.
 - Never obstruct the air inlets positioned on both sides of the device.
 - Always operate the device on an unencumbered hard surface.
 - Before every utilisation ensure the ventilation slits are not obstructed by any material.
 - Do not insert any objects into the ventilation slits.
 - Always unplug the power cable immediately after use.
 - The device is equipped with an easy to inspect safety fuse in case of malfunction. Before doing so, unplug the power cable. Repairs must be done only by authorised FLAEM personnel. Unauthorised repairs may be hazardous for the user.
 - Repairs must be done only by authorised FLAEM personnel. Any unauthorised repairs shall void the guarantee.
 - The average operating life of the compressor is 2,000 hours approximately.
- ATTENTION: Do not modify this device without authorisation from the manufacturer.**
- The Manufacturer, the Vendor and the Importer shall be held responsible for safety, reliability and performance only if: a) the device is used in compliance with the instructions for use b) the wiring where the device is being used is in compliance with safety regulations and current laws.
 - Interactions: the materials used which are in contact with medications have been tested with a wide range of pharmaceuticals. However, in view of the variety and continuous evolution of pharmaceuticals, interactions cannot be ruled out. We recommend the medication is used as soon as possible after opening and its prolonged exposure to the nebuliser is avoided.
 - The Manufacturer must be contacted to communicate unexpected problems/events related to functioning.

OPERATING INSTRUCTIONS

Before each utilisation, wash your hands accurately and clean the device as described in the paragraph “CLEANING, SANITISATION AND DISINFECTION”. During application it is advisable to suitably protect oneself from any dripping. The nebuliser and accessories are for personal use only in order to prevent any risks of infection.

It is recommended to replace them every six months.

1. Insert power supply cord (A8) into socket (A7) and then connect it to a mains outlet adequate to the device voltage. The position of the socket must be such that the device can be easily unplugged from the mains network

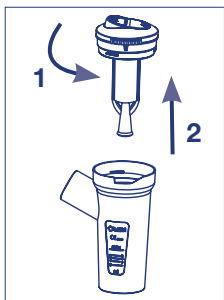
2. Open the nebuliser by turning the top (C5.1) anti-clockwise.



3. Pour the medication prescribed by the doctor into the bottom part (C5.3)

4. Close the nebuliser by turning the top (C5.1) clockwise.

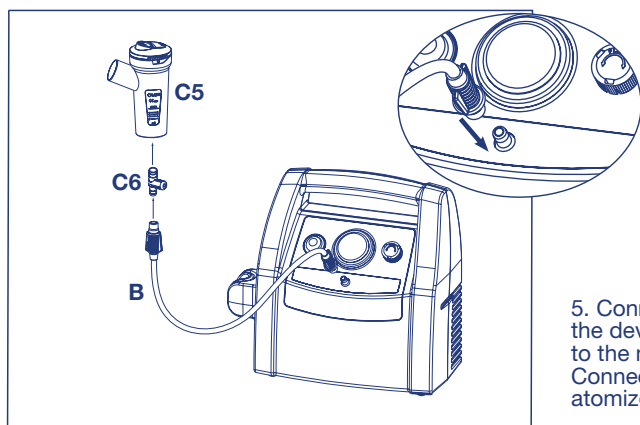5. Connect one end of the connection pipe (B) to the device. Connect the other end of the pipe (B) to the manual nebulisation control (C6). Connect manual nebulisation control (C6) to the atomizer (C5).

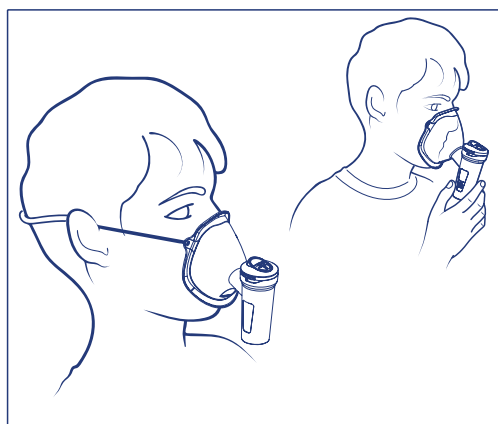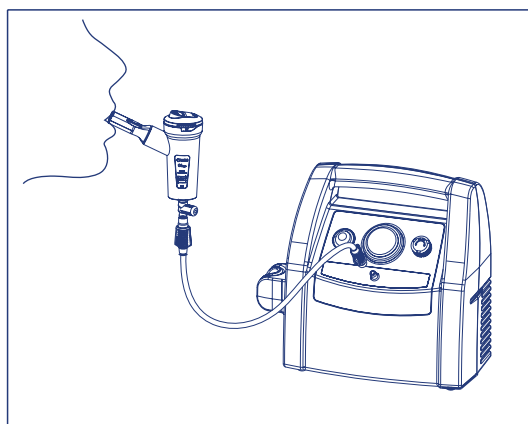6. Connect accessories as indicated in the "Connection diagram" on the cover. Sit comfortably holding the nebuliser in your hand, place the mouthpiece onto your mouth or alternatively use the nose piece or mask. Should you opt for the mask accessory, place it on your face as shown in the picture (with or without using the elasticated strap).

PRESSURE ADJUSTMENT


7.a Operating pressure



After connecting nebuliser RF7 Plus to the device with the Selector Button (C.4) closed in the position With Valve, operating pressure can be adjusted by rotating Knob (A3), clockwise to increase it or anti-clockwise to decrease it, until selecting the colour sector visible on the Pressure Gauge (A6) matching the desired particle size shown in the table below.

Throughout the device life, from time to time it is necessary to check maximum pressure attainable with the RF7 Plus nebuliser connected and with Knob (A3) turned on MAX.

If the maximum value of the green sector in the Pressure Gauge is not reached (A6), the device must be serviced.

RF7 PLUS TECHNICAL NEBULISATION FEATURES

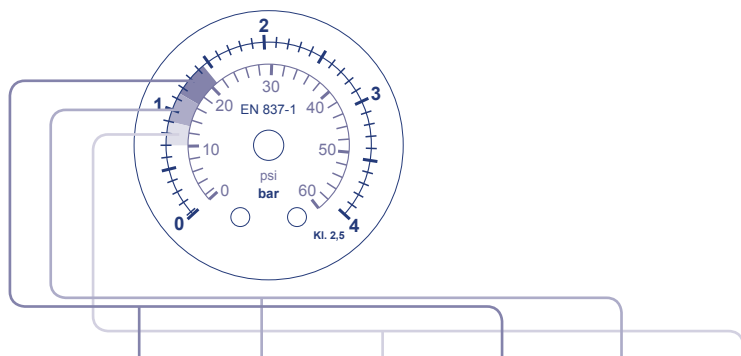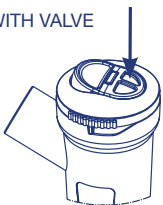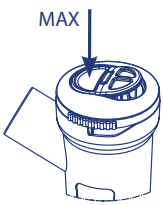	Green Area	Yellow Area	Orange Area	Green Area	Yellow Area	Orange Area
Pressure gauge value bar	1.15 ÷ 1.45	0.90 ÷ 1.15	0.70 ÷ 0.90	1.15 ÷ 1.45	0.90 ÷ 1.15	0.70 ÷ 0.90
Nozzle flow (1) l/min'	8.2	7.1	6.8	8.2	7.1	6.8
Selector Button Position	with Valve	with Valve	with Valve	MAX	MAX	MAX
MMAD µm (2)	2.21	2.53	2.44	2.38	2.72	2.95
(2) Breathable fraction < 5 µm (2)	84.2%	79.6%	80.7%	81.5%	77.8%	74.7%
Delivery ml/min' (1) approx	0.29	0.23	0.18	0.65	0.53	0.42
(1) Data recorded according to internal Flaem procedure I29-P07.5. Delivery Speed values are subject to variations depending on the patient's respiratory ability. (2) In vitro characterization certified by TÜV Rheinland LGA Products GmbH - Germany in compliance with European Standard for Aerosol therapy apparatuses EN 13544-1. Further details are available on demand.	WITH VALVE 			MAX 		

CHOICE OF DELIVERY METHOD

7.b Delivery method


Nebulisation can be obtained by maintaining the valve system active, or selecting the delivery mode at top speed, by means of the Selector Button (C5.4) of RF7 Plus atomiser. The mode must be selected after setting the operating pressure.

8. Start the device by means of switch (A1) and breathe deeply in and out.


9. Upon completing application, turn off the device and unplug it.

10. If after the therapy session an evident deposit of moisture forms within the pipe (B), detach the pipe from the nebuliser and dry it with the very ventilation from the compressor; this operation prevents possible blooms of mould inside the pipe.

SOFT TOUCH MASKS

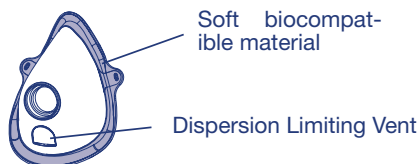
The new SoftTouch Masks feature an adaptive edge made of a soft biocompatible material ensuring an excellent adherence to the face and a Dispersion Limiting Vent. These innovations, which make them stand out from the rest, allow increased sedimentation of medication in the lungs and reduced dispersion into the air.

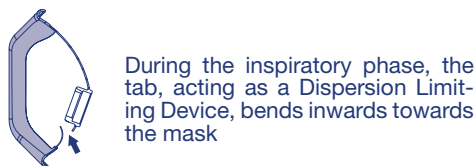
USE OF NEBULIZATION MANUAL CONTROL

To achieve continuous nebulisation action you should not use the manual nebulisation control (C6), especially in the case of children or persons with reduced physical, sensory, or mental capabilities. The manual nebulisation control is useful for limiting dispersion of the medication in the surrounding environment.

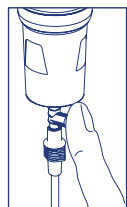To start nebulizing close with a finger the hole of the nebulizer manual control (C6) and breathe in gently; we recommend to hold your breath for a moment so that the inhaled aerosol droplets can be deposited,

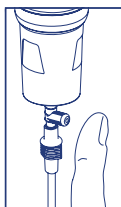meanwhile, to disable nebulizing, remove your finger from the hole of the nebulizer manual control to avoid the waste of drug, optimizing its acquisition. Then exhale slowly.

CLEANING, SANITIZATION AND DISINFECTION

Turn off the unit before any cleaning and disconnect the power cord from the mains.

UNIT AND OUTER SURFACE OF THE TUBE CLEANING

Only use a cloth dampened with antibacterial detergent (non-abrasive and free of solvents of any kind).

NEBULIZER AND ACCESSORIES

We recommend a personal use of both the nebulizer and accessories to avoid any risk of infection or contagion.

Open the nebulizer by rotating the upper part (C5.1) counter-clockwise, detach the nozzle by pressing as indicated by the 2 arrows in the "Assembly diagram" section C5.

SANITIZATION

Before and after each use, sanitize the nebulizer and accessories according to one of the methods describe below.

(Method A): The accessories that can be sanitized with the method A are (C1-C2-C3-C4-C5.1-C5.2-C5.3-C5.4-C6). Sanitize the accessories in warm (approximately 40° C) drinking water with mild dishwashing liquid (non abrasive) or in the dishwasher with hot cycle

(Method B): The accessories that can be sanitized with the method B are (C1-C2-C3-C4-C5.1-C5.2-C5.3-C5.4-C6) Sanitize the accessories by immersion in a solution consisting of 60% water and 40% white vinegar. Finally rinse thoroughly with warm (approximately 40° C) drinking water. Nach dem Reinigen schütteln Sie das After sanitization shake the accessories vigorously and place them on a paper towel, or alternatively, dry them with a jet of hot air (e.g. using a hair dryer).

DISINFECTION

The accessories that can be disinfected are (C1-C2-C3-C4-C5.1- C5.2-C5.3-C5.4-C6)


The disinfection procedure described in this paragraph is to be performed before using the accessories and is effective only if all the prescriptions herein are met and all the components involved are previously sanitized. The disinfectant to be used must be of electrolytic oxidizing chlorine type, (active ingredient: sodium hypochlorite), specific for disinfection and is available in all pharmacies.

Procedure:

- Fill a container, suitable to contain all the individual components to be disinfected, with a solution of drinking water and disinfectant, while respecting the proportions indicated on the packaging of the disinfectant itself.
- Completely immerse each individual component in the solution, taking care to avoid the formation of air bubbles in contact with the components. Leave the components immersed for the period of time indicated on the packaging of the disinfectant, and associated with the concentration chosen to prepare the solution.
- Remove the components now disinfected and rinse thoroughly with lukewarm drinking water.
- After disinfection shake the accessories vigorously and place them on a paper towel, or alternatively, dry them with a jet of hot air (e.g. using a hair dryer).
- Dispose of the solution as recommended by the disinfectant manufacturer.

AIR FILTERING

The device is fitted with a suction filter, located inside the pressure-adjustment knob, which must be replaced when it is dirty or discoloured. To replace it, lift the pressure-adjustment Knob (A3) and extract the Filter. Do not wash or reuse the same Filter. Regular filter replacement is necessary to assure proper compressor performance. The Filter must be regularly checked. Contact your retailer or authorised service centre for replacement filters.




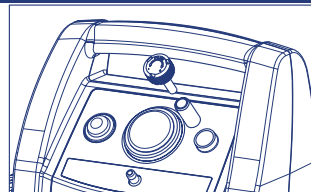
TROUBLE-SHOOTING

Switch off the device before any procedure and unplug the power cable from the socket.

PROBLEM	CAUSE	SOLUTION
The device does not work	The power cable has not been correctly inserted into the socket of the device or the power socket	Correctly insert the power cable in sockets.
The device does not nebulise or nebulises insufficiently.	The medication has not been inserted in the nebuliser.	Pour the right amount of medication in the nebuliser.
	The nebuliser has not been properly fitted.	Disassemble the nebuliser and reassemble it correctly as shown in the connection diagram on the cover.
	The nebuliser nozzle is clogged.	Disassemble the nebuliser, remove the nozzle and perform cleaning operations Failure to remove medication deposits from the nebuliser affects its efficiency and operation Strictly comply with the instructions contained in the CLEANING, SANITISATION AND DISINFECTION chapter
	The air pipe is not correctly connected to the device	Ensure the device's accessories are properly connected to the air intake (see connection diagram on the cover).
	The air pipe is bent, damaged or twisted	Unwind the pipe and ensure it is not crushed or punctured. If necessary replace it.
	The air filter is dirty	Replace the filter
The apparatus is noisier than usual	The pressure regulator is not correctly inserted in its housing	Fully and properly insert the pressure regulator in its housing

If after verifying the above mentioned conditions the device should not operate properly, contact your trusted retailer or the nearest authorised FLAEM service centre.

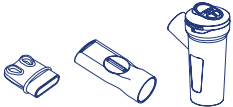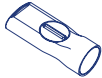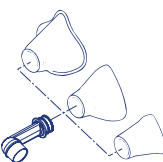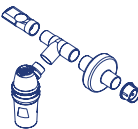SPARE PARTS AND ACCESSORIES	DESCRIPTION	CODE	PACKING
	RF7 Dual Speed plus nebulizer + mouthpiece with valve + nose piece	ACO464P	(1PC)
	RF7 dual speed plus nebulizer	ACO469P	(50PCS)
	2m pipe with fittings	ACO225P	(1PC)
	2m pipe with fittings	ACO457P	(50PCS)
	Mouthpiece with valve	ACO441P	(50PCS)
	Non-invasive nose piece	ACO448P	(50PCS)
	Adult mask	ACO473P	(50PCS)
	Paediatric mask	ACO472P	(50PCS)
	Complete newborn accessory set	ACO392P	(1PC)
	Newborn accessory set size 2	ACO391P	(1PC)
	Newborn accessory set size 1	ACO390P	(1PC)
	Newborn accessory set size 0	ACO389P	(1PC)
	PP adult nose piece	ACO450P	(50PCS)
	PP paediatric nose piece	ACO451P	(50PCS)
	Porex filter	ACO45	(50PCS)
	RF6P nebulizer kit w/ anti-dispersal filter	ACO388P	(1PC)
	Single-use filter for anti-leakage kit	ACO396P	(50PCS)

TECHNICAL DATA

Mod. P0712EM F3000 Plus Version

Power source:	230 V ~ 50 Hz / 210 VA
Fuse:	1 x T2A 250V
Max Pressure:	3,5±0,5 bar
Compressor air delivery:	14 l/min approx.
Noise (at 1 m):	55 dB (A) approx.
Operation:	Continuous
Operating conditions:	Temperature min. 5 °C; max. 35 °C Air humidity min. 10%; max. 95%
Storage conditions:	Temperature min. -25 °C; max. 70 °C Air humidity min. 10%; max. 95%
Operating/storage atmospheric pressure:	min. 690 hPa; max. 1060 hPa
Compliant with Directive 93/42 EEC	 0051
Apparatus size:	26 (W) x 12 (D) x 23,5 (H) cm
Weight:	2,4 kg

APPLIED PARTS

The applied parts of type BF are: (C1,C2,C3,C4)

TECHNICAL SPECIFICATIONS RF7 Dual Speed Plus Nebulizer

Minimum fill volume:	2 ml
Maximum fill volume:	8 ml

SYMBOLS

 Class II device

 Important: check the operating instructions

 Keep dry

I "ON" When the device is turned off, the switch interrupts compressor operation only on one of the two power supply phases.

O "OFF"

 Complies with: European standard EN 10993-1 "Biological Evaluation of Medical Devices" and European Directive 93/42/EEC "Medical Devices"

SN Serial number of device

 Type BF applied part

~ Alternate current

 CE Medical Marking ref. Dir. 93/42 EEC and subsequent amendments

 Hazard: electrocution.
Consequence: Death.
Do not use device while taking a bath or shower

 Manufacturer

DEVICE DISPOSAL

In conformity with Directive 2002/96/EC, the symbol shown on the device to be disposed of indicates that it is considered as waste and is therefore subject to "sorted waste collection". The user must therefore deliver (or have delivered) the above waste to an appropriate waste collection centre provided by the local authorities, or deliver it to the dealer when purchasing a new appliance of the same type. Pre-sorted waste collection and the subsequent treatment, recovery and disposal operations favor the production of appliances made of recycled materials and limit the negative effects to the environment and public health due to incorrect waste management. The unlawful disposal of the product by the user could result in administrative fines as provided by the laws transposing Directive 2002/96/EC of the European member state or of the country in which the product is disposed of.

ELECTROMAGNETIC COMPATIBILITY

This device has been designed to satisfy requirements currently required for electromagnetic compatibility (EN 60 601-1-2:2007). Electrical medical devices require special care. During installation and use with respect to EMC requirements, it therefore required that they be installed and/or used according to the manufacturer's specification. Potential risk of electromagnetic interference with other devices, in particular with other devices for analysis and treatment. Radio and mobile telecommunications devices or portable RF (mobile phones or wireless connections) may interfere with the operation of electrical medical devices. For further information visit www.flaemnuova.it or www.flaem.it Flaem reserves the right to make technical and functional changes to the product without notice.

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© Copyright by FLAEM NUOVA S.p.A., 2011 cod. 15612-0 rev. 06/2013